

BIOPAC Developer Suites



Scripting & Workflow • Network Data Transfer • Hardware API • Software API


BIOPAC Basic Scripting and Workflow

Scripting significantly reduces analysis time and improves consistency with standardized automated procedures. Workflow enables non-programmers to easily utilize scripting functionality.


Labs repeat the same procedure many times across multiple subjects. With scripting, it is possible to consolidate many tasks into one automated routine that eliminates the potential for error. This provides greater standardization and consistency while improving the quality of data.

Workflow's drag-and-drop interface allows researchers to automate and repeat the signal-conditioning and processing steps in an easy-to-use process without requiring programming or scripting knowledge. Workflow is included in all Scripting packages.

Scripting Tools

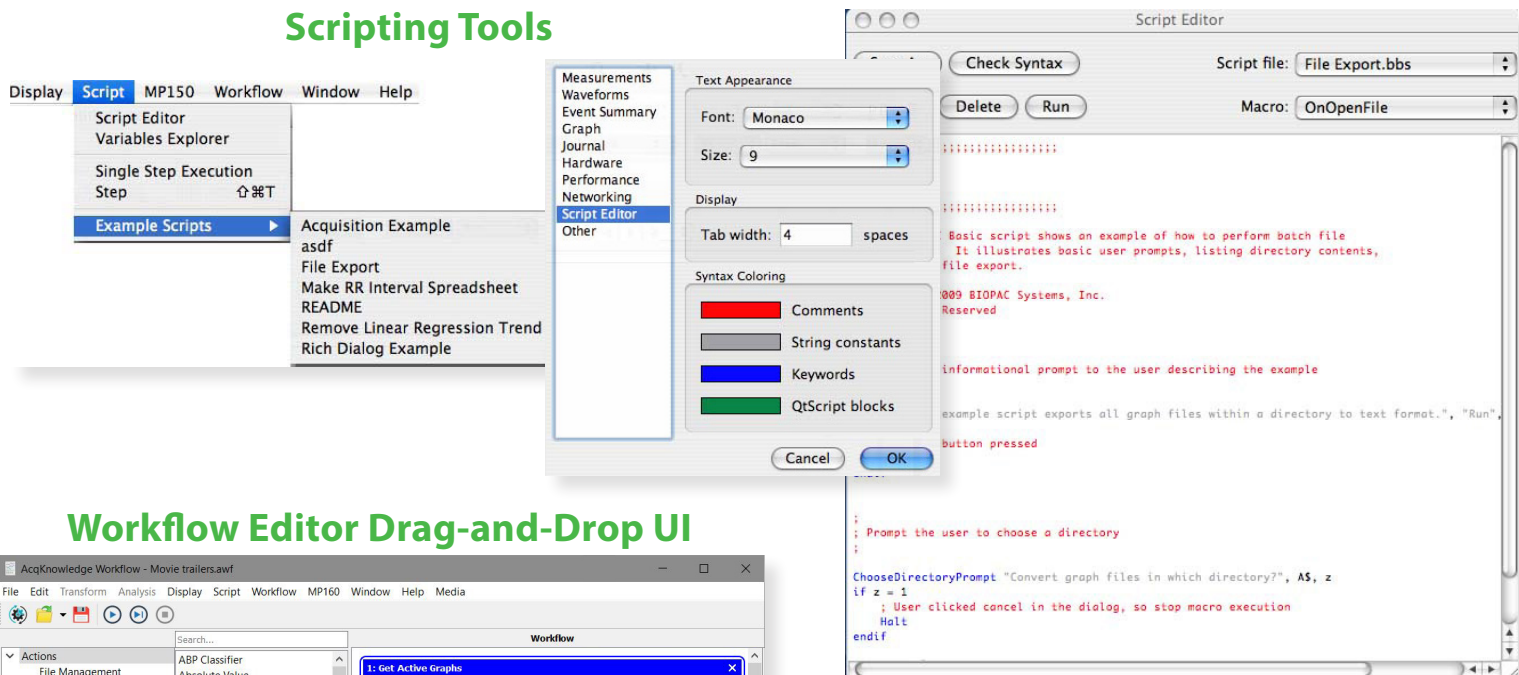

Workflow Editor Drag-and-Drop UI

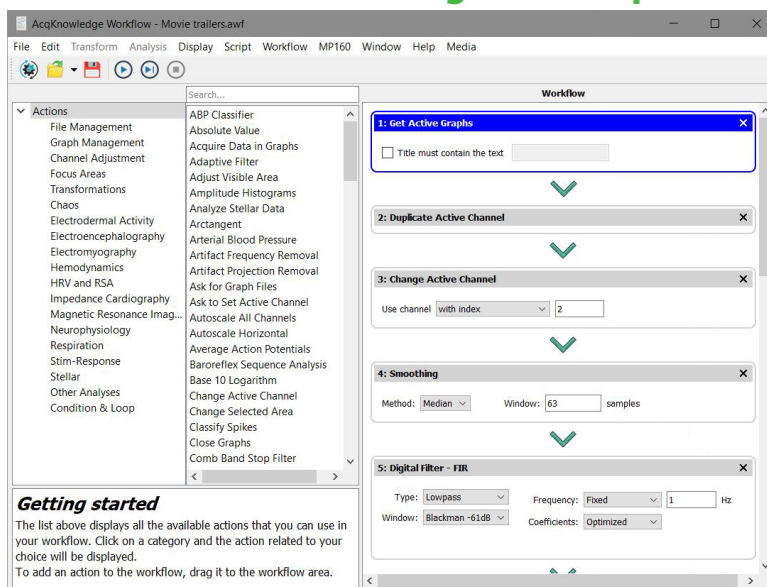

Customize acquisition setup and analysis

Automate acquisitions and perform batched analysis

Add dialogs and prompts to assist users

Improve quality with greater consistency

Expedite processing, reduce analysis time, and publish faster

BIOPAC

(805) 685-0066

info@biopac.com

www.biopac.com


Contact BIOPAC to learn more or request a quotation!

BIOPAC Developer Suites

Network Data Transfer

Network Data Transfer (NDT) is a real-time data transfer system that allows access to the data being acquired into a graph by *AcqKnowledge* for use in an external application. The computer used to run the *AcqKnowledge* process and the custom application may be the same computer.

- Allows third-party applications—including the Vizard VR development environment—to tap into the data stream generated by the MP unit and *AcqKnowledge* during acquisitions.
- Provides networking facilities that allow for integration into a distributed application environment.
- Provides basic controls to allow external applications to query and control the *AcqKnowledge* application state.

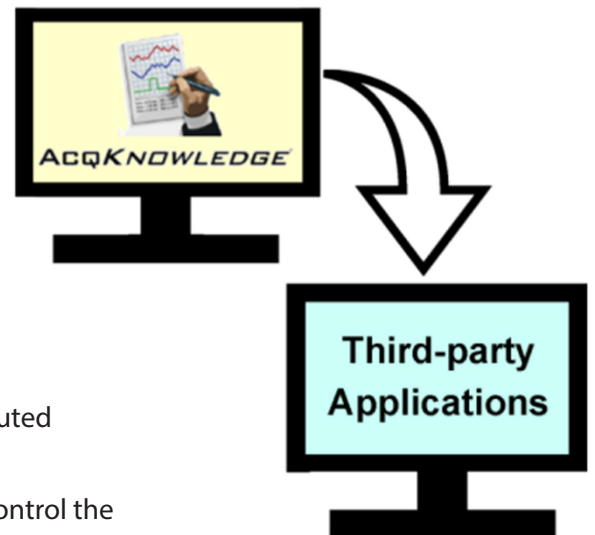

Hardware API Research or Education MP Systems (Windows only) **Not available for MP40 or MP45**

The hardware API gives developers control over BIOPAC MP data acquisition units to do the following:

- Acquire data
- Acquire at different sample rates
- Set triggers
- Get status of the MP unit
- Use Analog Output channels
- Use Digital I/O

The implementation of these functions is compiled into a Windows 32-bit DLL called "mpdev.dll." The interface is documented in C/C++, but any programming language that can utilize Windows 32-bit DLLs can access the BIOPAC Hardware API. The hardware API includes sample projects for C/C++, C#, LabVIEW (v8.2), MATLAB, and VB.NET.

Site licenses are available for the Hardware API; see BHAPI-9, BHAPI-24, and BHAPI-99 to extend the number of concurrent licenses.

Software API read-only access to *AcqKnowledge* 4.4 or above generated files (Windows only)

The BIOPAC File Format Application Programming Interface (API) is a software library to identify and parse information in BIOPAC's ACQ binary file format for programmers to use when creating applications for alternative analysis. Use the base functions in a variety of combinations for a new application.

- Initialize an ACQ file structure
- Close an ACQ file structure
- Retrieve channel information
- Retrieve samples by segment of a specified channel
- Retrieve all samples of a specified channel
- Retrieve a particular sample of a specified channel
- Retrieve samples by time slice of a specified channel
- Retrieve Journal text
- Retrieve marker information
- Retrieve marker text of a specified marker



Video Tutorials Available!

Contact BIOPAC to learn more or request a quotation!

BIOPAC

(805) 685-0066

info@biopac.com

www.biopac.com